

# Invacare® Matrix® | Mini Back

## Features

**Cover** – Foam and cover overlap the shell for comfort and protection. Water resistant and breathable, available in a choice of neon pink, bright blue or Kid\*ab\*ra print.

**Hardware** – “Mini Set” hardware installs and adjusts in minutes. Small attachment point for minimal interference on back canes. Easily adjusted for growth – depth (4”), angle (40°) and even rotation. Adjustable mounting pins allow fit to chairs same or 1” wider than back. Quick release for convenient portability.


**Back shell** – Light weight, durable aluminum shell, designed for easy attachment of optional lateral supports and headrest.

**Adjustable Lumbo-Sacral Support** - Included.

**Headrest Adapter** - The back shell is per drilled to accept Motion Concepts, Stealth and Otto Bock™ headrest mounts.

**Swing-Away Lateral Supports** - 2.5” of vertical adjustment on back pan. Mounts inside back shell without interfering w/depth/angle adjustment Easy push button release, available in Black, Blue & Pink.

**Lateral pads** - XS: 3.25” H x 4.25” L

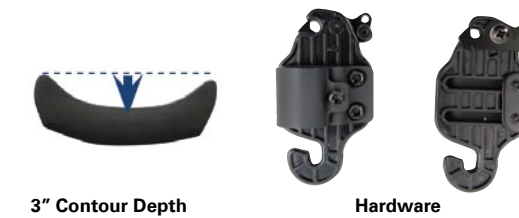

3” Contour Depth

Hardware

## Technical Data

Height	Width	Depth Adjustment	Angle Adjustment	Product Weight	Weight Capacity	Limited Warranty
7", 10", 12", 14" 18, 25, 30.5, 35.5cm	10" - 14" 25 - 35.5cm	4" 10cm	40° range +/- 20°	1.3lb and up 0.5kg and up	150lb 68kg,	2 years



### Designed with kids in mind

The **Invacare® Matrix® Mini Back** features Mini Set hardware, a miniature version of the popular Easy Set hardware, designed with kids in mind. The Mini back is contoured to fit kids of all shapes and sizes – combining exceptional support in an appealing low profile package.

**Codes:** Back - MINlxxxx  
Cover - MCxxxx

