

PS Cushion

SUPERIOR POSITIONING, STABILITY AND COMFORT



The Invacare® Matrix® PS Cushion design features an anatomically contoured shape with a unique waffled ischial relief recess and an innovative reversible outer cover. It is a great choice for a very wide range of individual needs, including pediatric, geriatric, bariatric and highly active.



matrix[®]
SEATING SERIES

PS Cushion



Outer cover:

- Innovative reversible Startex™ outer cover may be used smooth or fabric side out. Hook and loop fastener strips on both sides.
- Startex™ fabric is moisture resistant and breathable

Inner cover:

- Moisture resistant
- Easily wiped clean

Cushion:

- Durable high resilient foam prevents “bottoming out”
- Anatomically molded shape for optimal positioning and comfort
- Waterfall front edge
- Made with Ultra-Fresh™ for antibacterial and odor protection
- Waffled recess provides immersion of body prominences, helps maintain pelvic position and prevents sliding

Options




Cushion and Cover Modifications:

Standard modifications include seat width and depth reductions, as well as leg length discrepancies



Pelvic Obliquity Kit:



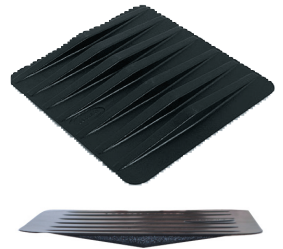



Designed to accommodate up to 1° obliquities.



Cushion Rigidisers:

Lightweight plastic seat insert designed to fit to the bottom of the cushion inside the cover. Sits between seat rails for level support and eliminates hammock effect of sling seat upholstery.

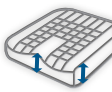Tech Specs

	 Width*	 Depth**	 Height	 Product weight	 Weight Capacity	 Limited Warranty
Paediatric	8-14"	10", 12", 14", 16"	3.75" Adductor 3.5" Abductor 2.5" Leg trough	From 450g (1lb)	68kg (150lb)	2 years
Standard	14, 16, 17, 18, 20"	16, 18, 20, 22"	4" Adductor 3.5" Abductor 2.5" Leg trough	From 900g (2lb)	136kg (300lb)	2 years
Heavy Duty	21" - 30"	18" - 30" 40.5 - 76cm	4"/3.5"	From 3kg (7lb)	272kg (600lb)	2 years

* Actual width is 1/2" narrower than stated width

** Actual depth is 1 1/2" longer than stated depth



FOR FURTHER INFORMATION:

New Zealand: Freephone 0800 468 222, email sales@invacare.co.nz or visit www.invacare.co.nz

Australia: Freephone 1800 460 460, email sales@invacare.com.au or visit www.invacare.com.au

Asia: Phone +61 2 8839 5330, email exportsales@invacare.com.au or visit www.invacare.com

Thailand: Phone +66 0 2108 8555 or visit www.invacare.com



Yes, you can.®